


OP-D-OP Polysoft Face Shield

COMPONENTS

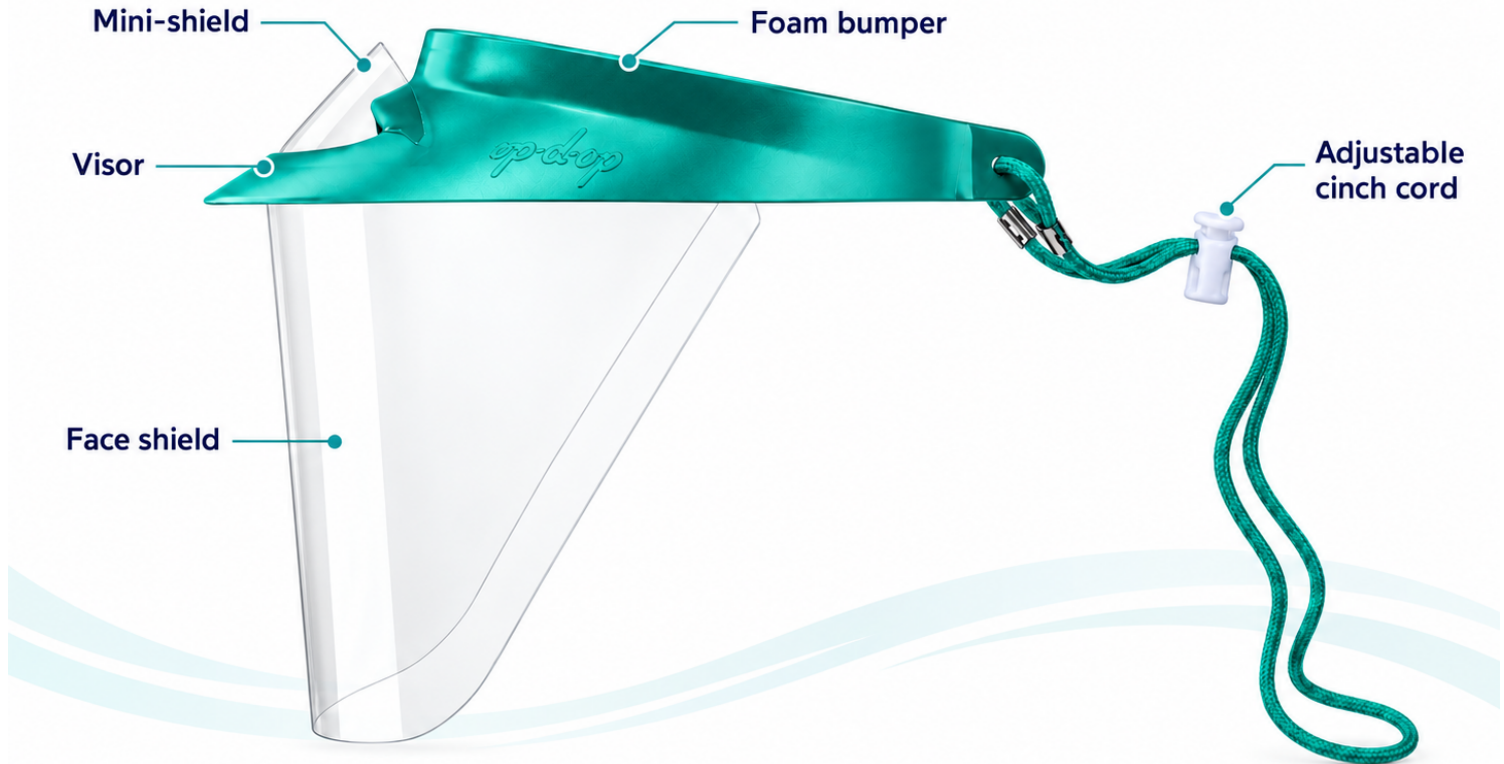

Supplied by Veliro Medical Supplies - Authorised OP-D-OP Reseller

IMPORTANT

Always wear your OP-D-OP with the face shield and mini-shield attached.

Replacement Shields and Foam

Replacement face shields, mini-shields and foam bumpers are available as needed. Visors should only be used with shields manufactured by OP-D-OP.

Use & Safety Information


OP-D-OP face shields are designed to be worn over the face of health care workers to reduce the risk of being splattered or splashed with possibly infectious bodily fluids. OP-D-OP products do not prevent the contraction of disease or infection and should always be worn using standard precautions.

Warning: Materials which may come into contact with the wearer's skin could cause allergic reactions to susceptible individuals.


Instructions for Use

OP-D-OP Polysoft Face Shield


1 REMOVE MINI-SHIELD

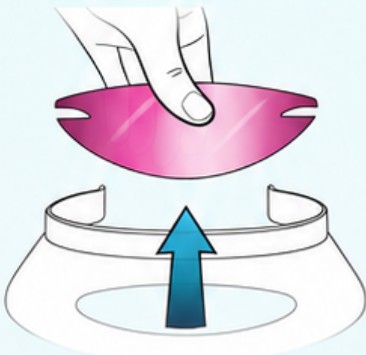


2 PEEL SHIELDS

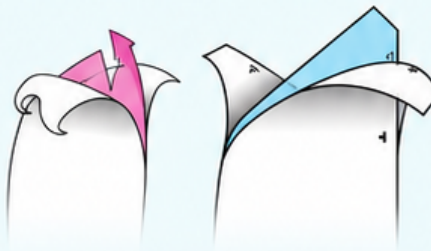

3 ATTACH FACE SHIELD

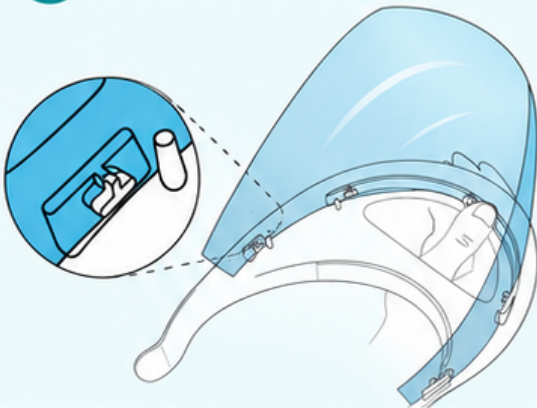

4 RE-ATTACH MINI-SHIELD

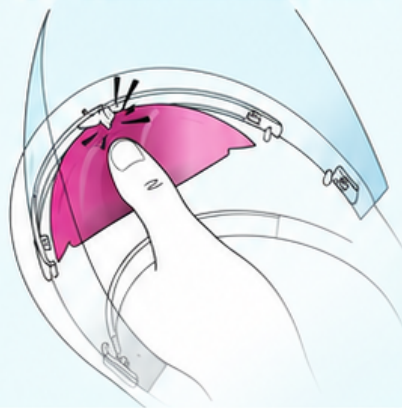

5 WEAR IT

To prevent fogging, wear your Op-d-op with the shield parallel to your face – fit two or three fingers between your chin and the shield to check.

You're ready to go!



Storage



Store and transport your Op-d-op in its original packaging until use.



Keep the protective film on the shield during transport.



Replace visors after 5 years from the date of purchase.

Cleaning your OP-D-OP



Warning: If your shield is heavily splattered or scratched, properly dispose of it and replace with a new one.



To reuse shields, wash thoroughly under running water using a suitable antimicrobial soap (e.g. Hibiclens or Bactoshield). Rinse clean and dry with a soft towel.



Shields can also be cleaned using a cold sterilization solution or alcohol (ethanol 50–80%). Always follow the manufacturer's sterilizing instructions.



Do not autoclave shields or visors.



Always keep the shield assembly right side up when cleaning.

